


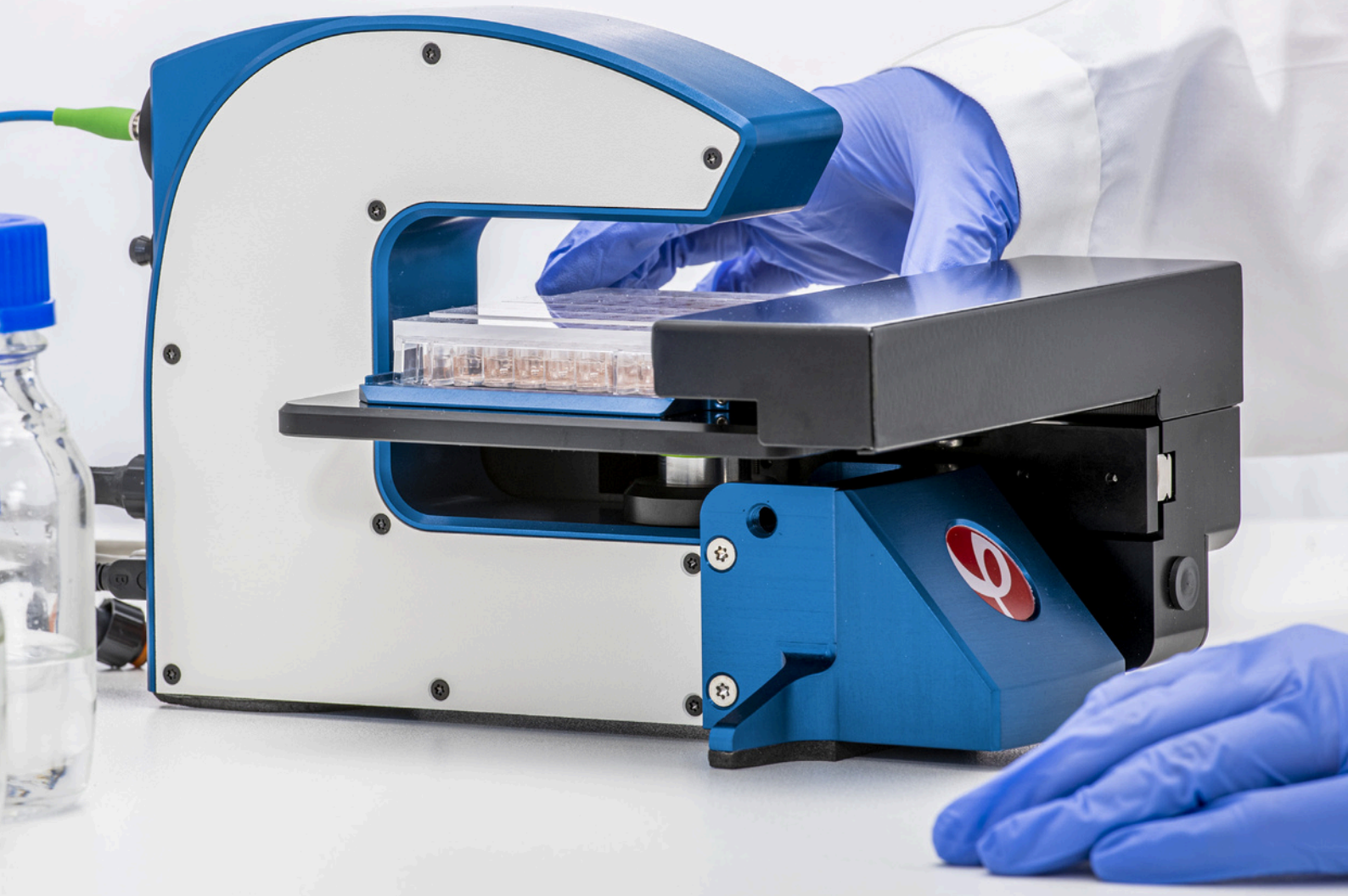


HoloMonitor®M4


Non-invasive Live Cell Analysis

The HoloMonitor® live cell imaging system enables long-term non-invasive analysis of cell cultures within your standard incubator. Cell biologists worldwide use label-free cell imager to automatically obtain live-cell morphology, migration, and proliferation data down to a single-cell level in real-time.


Focus on Your Cells

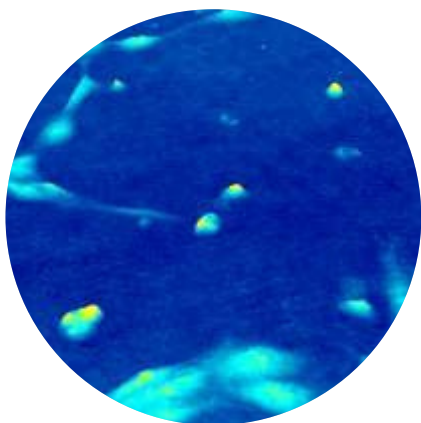

Induced pluripotent
stem cells (iPS)

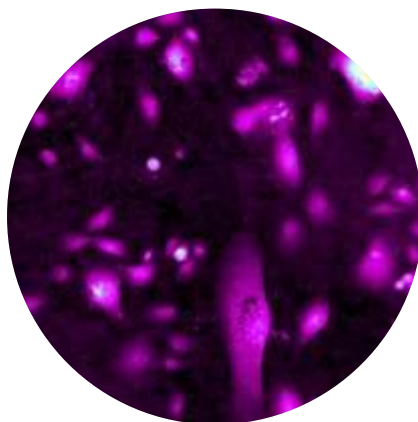


Primary human small airway
epithelium cells (SAEC)

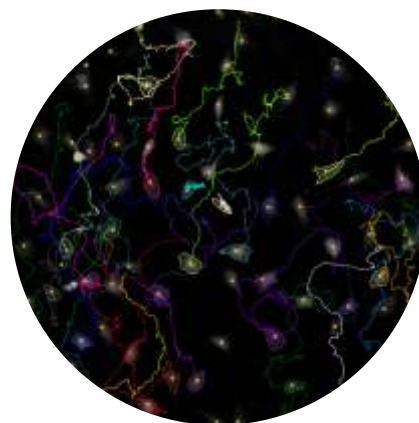


Single cell tracking

**Non-invasive
Quantitative
Single cell analysis
Fits inside the incubator**

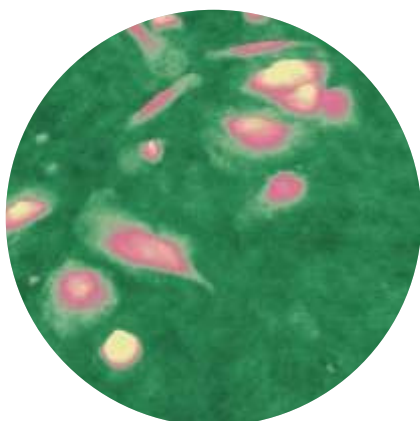

Macrophages

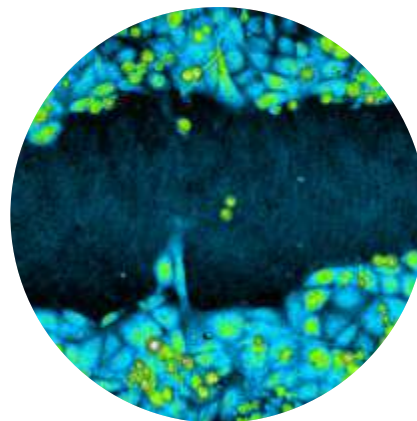

Wound healing of breast



HoloMonitor® is a live cell imaging system designed to visualize and quantify adherent cells over time. It is gentle to all cell types and can be used to monitor cell lines as well as delicate cell types, such as induced pluripotent stem cells or primary human cell cultures. As cell imaging with HoloMonitor is label-free and harmless to your cells, the cells can be used for further studies after imaging.

Our users have employed the system in a variety of fields, including cancer research, drug discovery or chemical toxicity studies. Research results are published in more than 200 peer-reviewed publications. All publications can be found on our website www.holomonitor.org


User-friendly Applications

The HoloMonitor software (App Suite) is a powerful tool designed to extract quantitative data about your cells from the holographic images captured by the instrument.

The App Suite software is intuitive and guides the user through the simple step-by-step work flow, from experiment setup to result analysis. You can also easily create colorful videos and images for peer-reviewed publications. All data can be exported to Excel for further analysis.


Guided End-point Assays

Ensure experiment reproducibility and optimal cell quality by checking up on your cells before setting up a time lapse experiment.

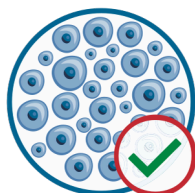

CELL QC ASSAY

Evaluate cell integrity and detect changes of your cell culture.

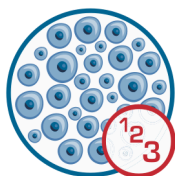

CELL COUNTER

Determine the exact cell number in your cell suspension.




Guided Kinetic Assays

Tailored kinetic assays that automatically analyze images of your cells and provide kinetic data on your cell populations behavior in real-time.




DOSE RESPONSE ASSAY

Recognize drug toxicity patterns and kinetics of drug effects from interactive dose response curves.

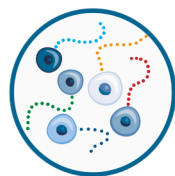


CELL MOTILITY ASSAY

Explore average cell population motility based on cell speed and accumulated cell distance in real-time.

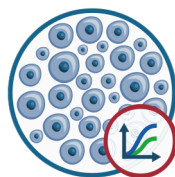

CELL PROLIFERATION ASSAY

Robust kinetic cell proliferation assessment by direct identification and counting of individual cells.

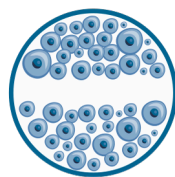

WOUND HEALING ASSAY

Easily obtain kinetic information on how cell populations close the wound gap and determine cell front velocity.


In-depth Analysis

The single-cell tracking functionality enables in-depth analysis of both individual cells and entire cell populations at the same time




SINGLE CELL TRACKING

Study single cell motility and migration patterns and distinguish between cell subpopulation behavior

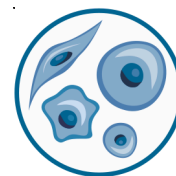

KINETIC MORPHOLOGY


Examine more than 30 parameters, as cell morphology can be sensitive markers for cell damage.

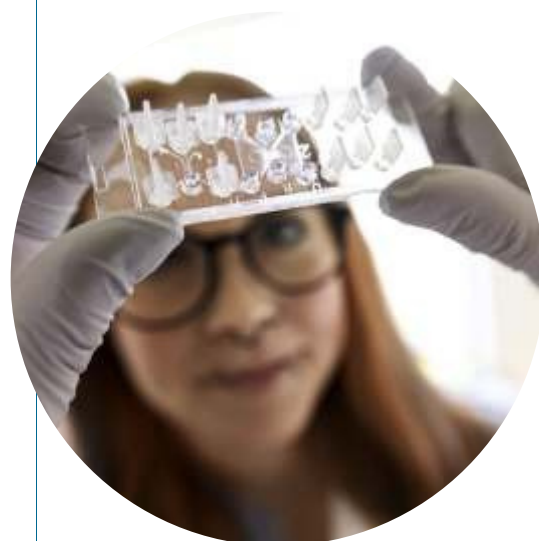

Standard Way vs. HoloMonitor® Way

Standard cell culture methods require various instruments, from plate readers to microscopes. Experiments require multiple setups and consumables and most readouts are indirect as label intensities are measured.

HoloMonitor®, on the other hand, enables real-time label free and quantitative studies of your cells directly in the incubator. All information is directly obtained by monitoring cells without any phototoxic damage or other disturbances.

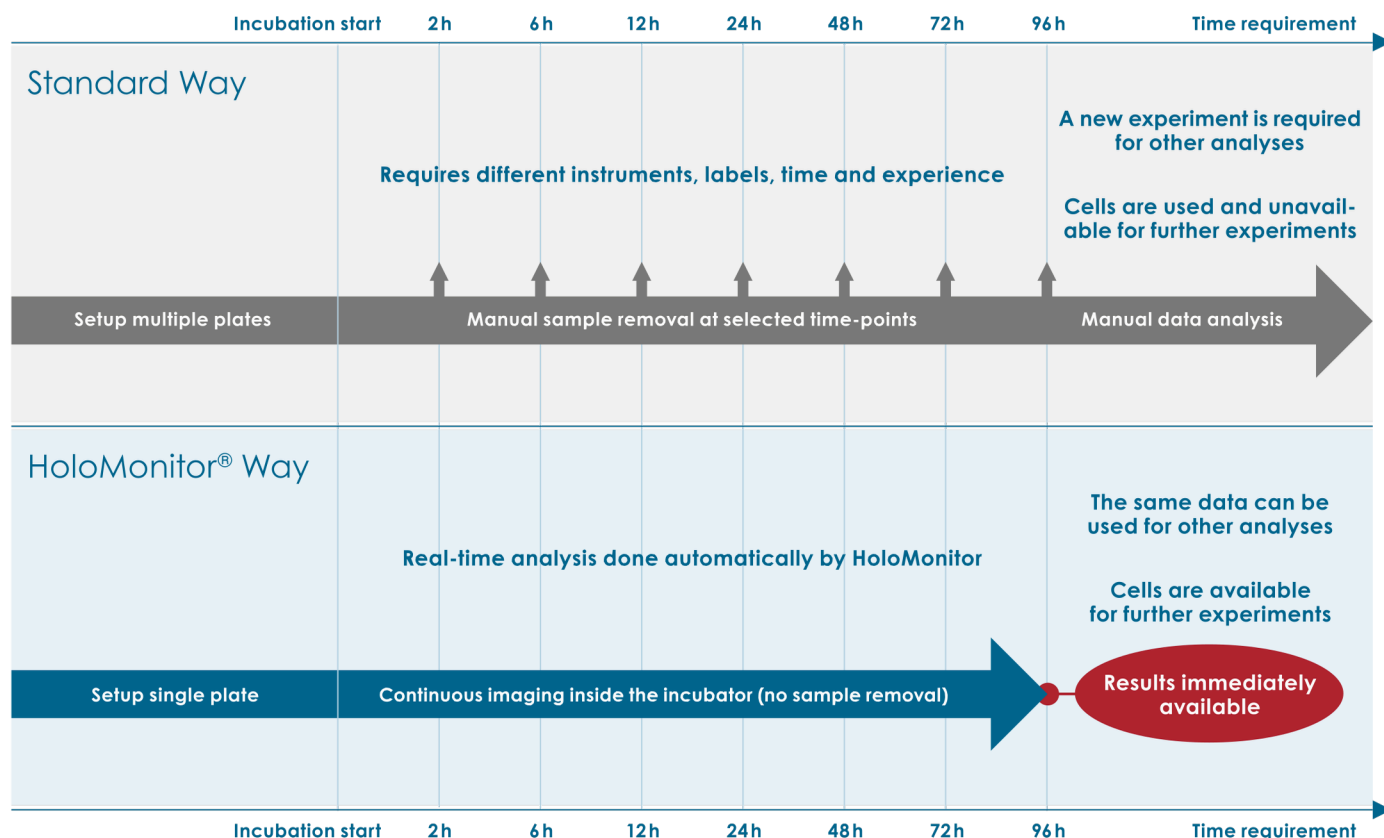

HoloMonitor®M4 Features


- Incubator tolerant design – cells are kept undisturbed in their preferred environment during the entire experiment.
- Harmless laser illumination power – enables long-term time-lapse image acquisition and preserves your cell integrity.
- Motorized stage ensures precise multiple position imaging – capturing the full story of your cells.
- Easy to use, powerful and unique software visualizes and provides quantitative data about your precious single cells or cell populations.
- Great image quality – carefully designed vessel lids, HoloLids™, for perfect images.

HoloMonitor®M4 Benefits

- Save hands-on lab time, money and cells
- No expensive consumables or vessels needed
- Label-free technology with no phototoxicity
- Minimize user bias with automated analysis
- Monitor your cells 24/7 and get results in real-time
- Get direct and quantitative measurements
- Re-use your imaged cells further for other assays
- Re-analyze data for multiple results from one experiment


Defining Your Cells Without Labels

HoloMonitor® is based on quantitative phase imaging enabling non-invasive visualization and quantification of living cells without compromising cell integrity. Unlike phase contrast microscopy, individual cells can be identified using the HoloMonitor technology taking your research to new discovery horizons.

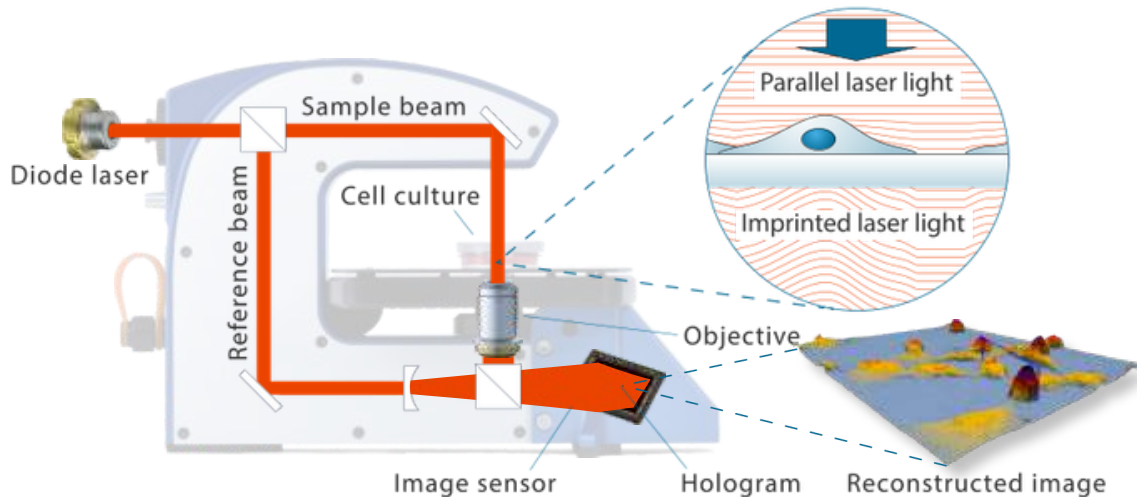


Holographic microscopy creates quantitative phase images by letting a sample beam and a reference beam interfere to generate an interference pattern – the hologram.

The hologram is recorded by an image sensor and processed by the computer to produce a reconstructed three dimensional (3D) image.

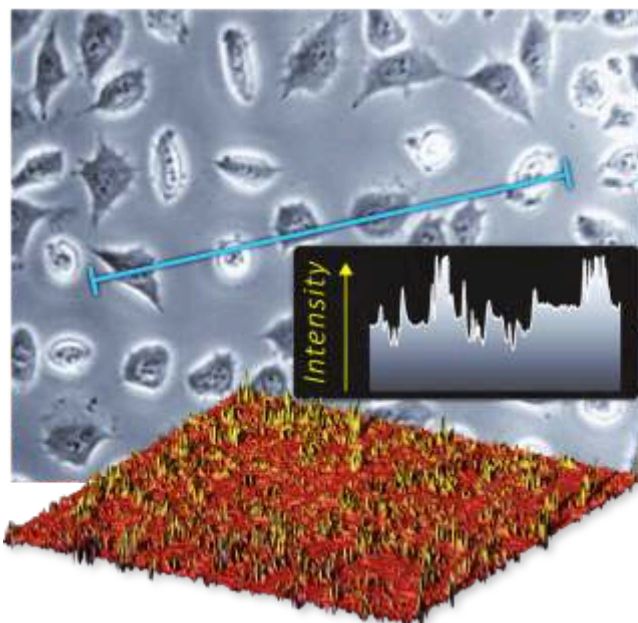

Phase Contrast Image

- Using phase contrast, image background intensity values cannot be accurately determined.
- It is impossible to segment cells or display intensity as 3D image.

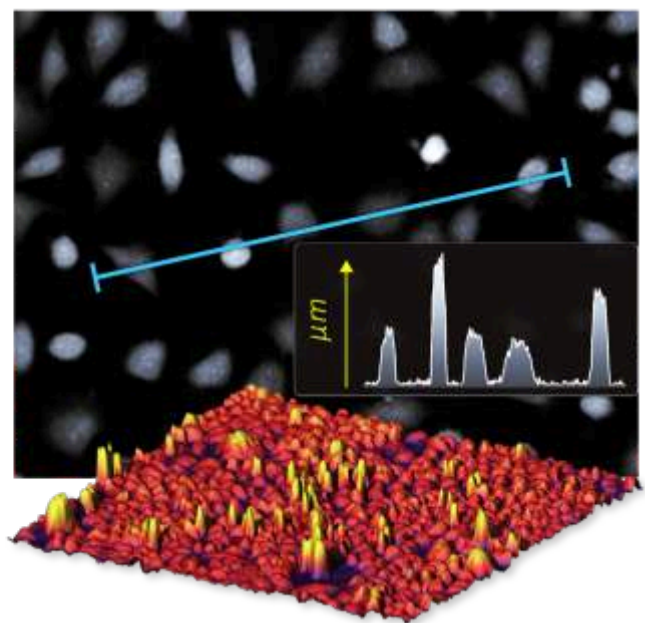

Holographic Image

- Holographic images have a background level of zero.
- It is easy to distinguish cells in the 3D quantitative phase image, as each single cell creates its own peak.